




POLYPRO PREFABS

Pre-stressed, semi rigid polypropylene prefabricated devices that incorporate many features of custom orthoses at a fraction of the price.

Polypro Prefabs are ideal as an economical first step in the treatment of a variety of common foot, ankle and lower leg pathologies that are secondary to mild to moderate subtalar joint pronation, excessive mid-tarsal joint motion and/or mild forefoot deformities.

Polypro Prefabs are also an excellent choice as an inexpensive "second pair" for those patients with limited insurance coverage.

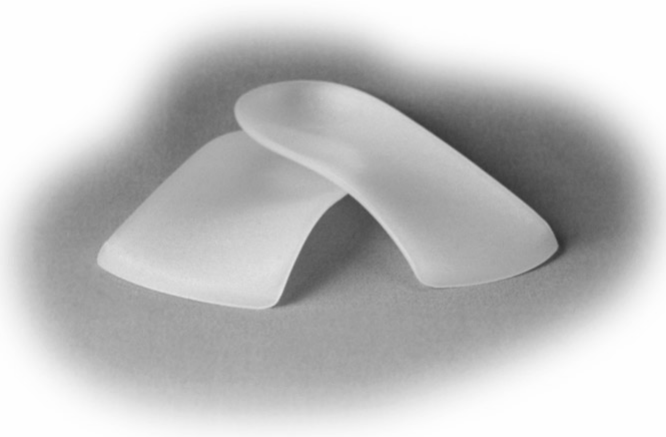

Polypro Prefabs were designed by ProLab Orthotics/USA, one of America's most respected prescription foot orthoses labs, and incorporate a number of advanced functional control features including:

- 2mm medial heel skive to enhance rear foot control and decelerate subtalar joint pronation,
- 3 degree intrinsic forefoot valgus correction to load the lateral column and limit mid-tarsal joint motion,
- Traditional "functional" shell shape parameters,
- 14mm heel cup, and
- A textured surface that eliminates the need for a top cover.

Polypro Prefabs are sized by shoe size and are available in men's and women's sizes 6 through 13. They are sold in single units as well as an 18 pair (mixed size) starter package.



PROLAB
ORTHOTICS/USA

For more information or to request an order form please call:

Paris orthotics

Paris Orthotics Ltd.
3630 East 1st Avenue
Vancouver BC V5M 1C3
Toll Free: 1.800.848.0838

Paris

WALKWELL 328 UNIVERSAL POST OPERATIVE SHOE



WalkWell's new 328 Universal Post Operative Shoe is a welcome alternative to the many poorly constructed, ill-fitting post-op shoes currently on the market.

The **WalkWell 328** Features

- Rigid rocker sole that effectively offloads metatarsal head pressure by eliminating metatarsal-phalangeal joint dorsiflexion,
- Bacteria static EVA insole that comfortably attenuates shock,
- Full reinforced heel counter with a unique adjustable closure system designed to accommodate heels of almost any size or shape,
- Soft and breathable micro-fabric upper with three Velcro closures to gently accommodate edema, post-op dressings or accommodative orthoses.
- Aesthetically pleasing all black design that promotes patient compliance.

Clinical Indications

- Osteotomies with and without fixation
- Post nail surgery
- Metatarsal fractures
- Metatarsal phalangeal joint injuries
- Metatarsal phalangeal synovitis and tendonitis
- Hallux limitus
- Hallux Valgus
- Neuropathic ulceration
- Persistent edema

The **WalkWell 328** rigid rocker sole design is the best in its class. A study published in the Journal of the American Podiatric Medical Association indicates that WalkWells can reduce peak forefoot pressures up to 20% better than conventional post-op shoes.

WalkWell 328's are universal (fit left or right) and are available in seven sizes from Women's Small (5) to Men's XX-Large (16.5). 328's are sold in single units as well as a 14 shoe (mixed size) starter package.

For more information or to request an order form please call:

Paris orthotics

Paris Orthotics Ltd.
3630 East 1st Avenue
Vancouver BC V5M 1C3
Toll Free: 1.800.848.0838

